

SAFE PERCUTANEOUS DRAINAGE KIT

For all-purposes drainage procedures

D R E N T E C H
UNICO  **TM**



STANDARD



MULTI




J SELDINGER



FORTY

PULL & PUSH

Thanks to double integrated valve drainage procedure can be completed by mean of simple "pull and push" manoeuvres. Three-way stopcock no longer needed.

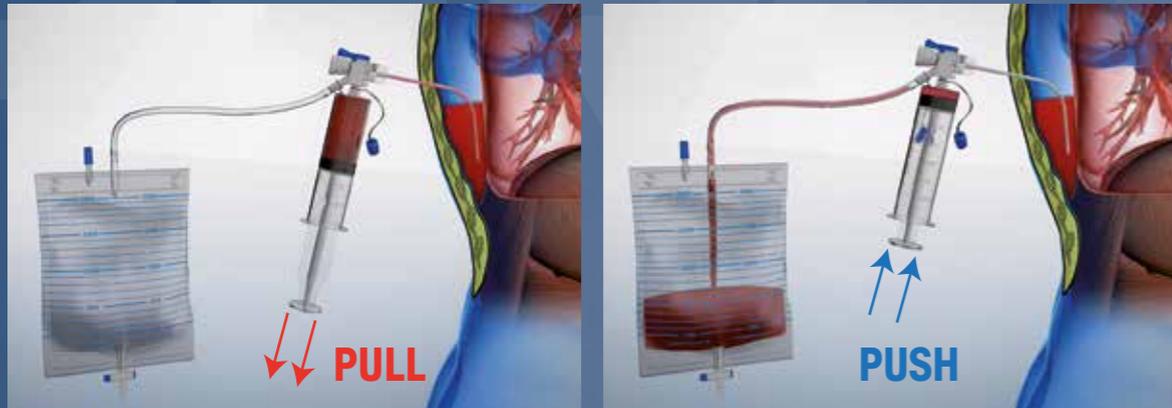

EXTENSION LINE

By connecting an extension line it improves patient comfort in case of extended use of the device and more flexibility for users.

VERRES NEEDLE

retractable tip to achieve safe and efficient access. The blunt tip retracts during penetration and extends when the needle enters free space.



POLYURETHANE CATHETER:

totally radiopaque, anti kinking and cut-resistant. It is biocompatible and thermosensitive which makes it the first choice for long-term drainage. Larger inner diameter increases drainage flow rate.

DOUBLE INTEGRATED ONE-WAY VALVE

designed to improve comfort and safety of your drainage procedures. The valves ensure patient and users safety by reducing risk of external environment contamination, pneumothorax and fluid leakage.

Veress needle with luer lock port for primary positioning check.

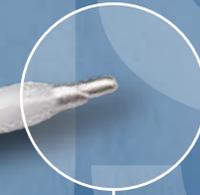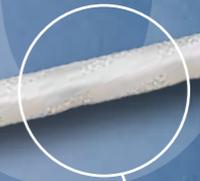

The best for pneumothorax treatment (primary and secondary), malignant or chronic pleural effusions and uncomplicated empyema.



STANDARD

	CODE	PCS/BOX
Percutaneous drainage Kit - CH/Fr 9	10800	6
Percutaneous drainage Kit - CH/Fr 12	10801	6
Percutaneous drainage Kit with extension line - CH/Fr 9	10802	6
Percutaneous drainage Kit with extension line - CH/Fr 12	10803	6
Percutaneous drainage Kit with extension line - CH/Fr 15	10830	6

EXTRA LONG*


	CODE	PCS/BOX
Percutaneous drainage Kit with extension line XL- CH/Fr 9	10820	6
Percutaneous drainage Kit with extension line XL- CH/Fr 12	10821	6

*50% Longer than standard

Particularly suitable for abscess, nephrostomy pericardiocentesis and miscellaneous drainage procedures.



PIG TAIL

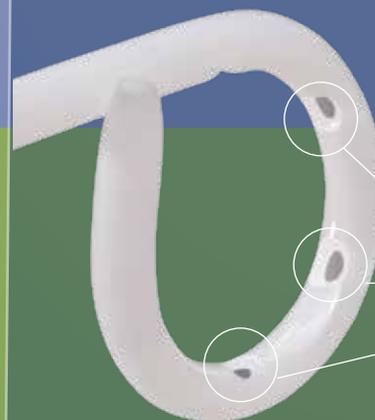

Pig tail provided with five inner drainage holes, to prevent irritation of surrounding tissues and occlusion of the catheter

MULTI


	CODE	PCS/BOX
Percutaneous drainage Kit "MULTI" CH/Fr 9	10816	6
Percutaneous drainage Kit "MULTI" CH/Fr 12	10817	6
Percutaneous drainage Kit "MULTI" with extension line - CH/Fr 9	10818	6
Percutaneous drainage Kit "MULTI" with extension line - CH/Fr 12	10819	6
Percutaneous drainage Kit "MULTI" with extension line - CH/Fr 15	10831	6

The most precise positioning through a fast-track Seldinger technique, all together with the unmatched UNICO's features.



J SELDINGER

	CODE	PCS/BOX
UNICO-J Percutaneous Drainage Seldinger - CH/Fr 9	10823	6
UNICO-J Percutaneous Drainage Seldinger - CH/Fr 12	10824	6

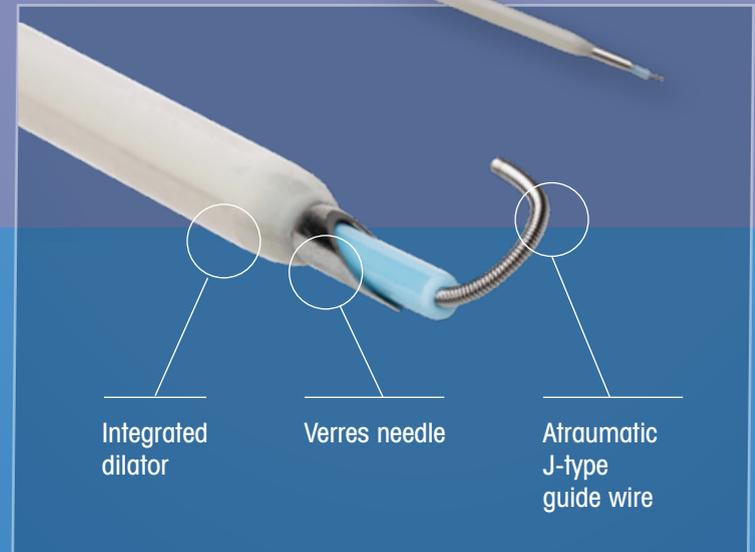

D R E N T E C H
UNICO™
FORTY

The forty centimeters percutaneous drainage design for apical positioning to ease liquids and air evacuation

FORTY

	CODE	PCS/BOX
Percutaneous drainage Kit with extension line - CH/Fr 9	10810	6
Percutaneous drainage Kit with extension line - CH/Fr 12	10812	6

D R E N T E C H
UNICO™
BASE


ONE SHOT PUNCTURE SET
 For pleural drainage of liquids and air, particularly suitable for: pneumothorax, effusion and low viscosity empyema

BASE


	CODE	PCS/BOX
Unico Base - gauge 13/ 2,5 mm	10815	20

ACCESSORIES

	CODE	PCS/BOX
Drainage bag	12019	25
Catheter fixing kit	10546	50
M connector for chest drain	10548	50
F connector for chest drain	10504	20

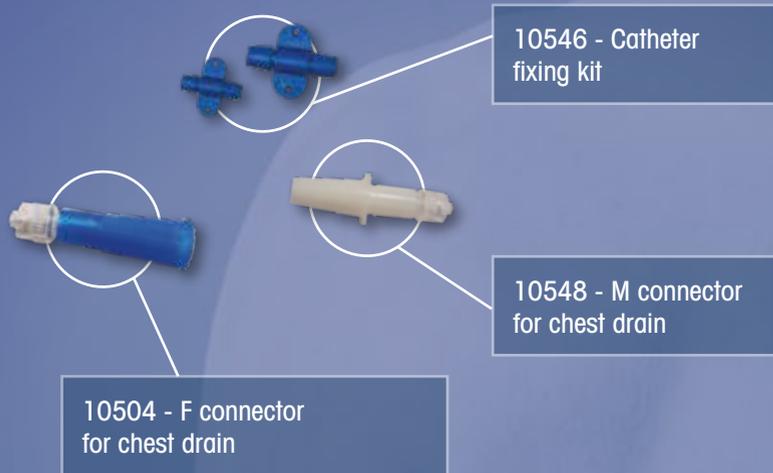

INTEGRATED ONE-WAY VALVE
to prevent liquid back flow

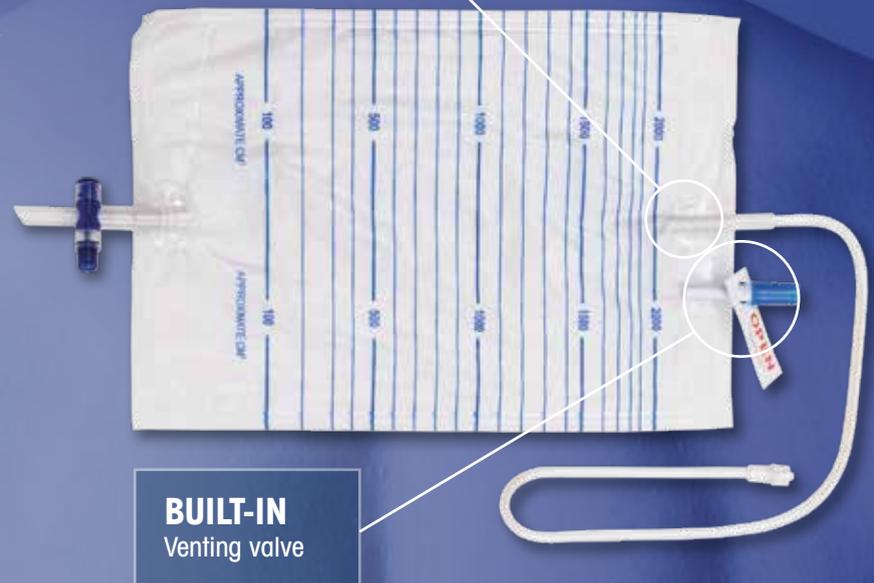

REDAX™

REDAX - Via Galileo Galilei, 18
46025 Poggio Rusco (MN) ITALY
Tel. +39 0386 830582 Fax +39 0386 51898
E-mail: info@redax.it - www.redax.it

UNICO
video tutorial

