

DRENTECH™ PALMEVO

Complete thoracic drainage system
for suction, air leakage monitoring and intrapleural data

SUCTION MODE

- Totally independent from any hospital vacuum sources

- Suction regulation
-10 up to -60cmH₂O (-1 to -6 kPa)

- Portable vacuum unit with rechargeable battery
up to 36 hours of use

MONITORING MODE

- Real-time and historical monitoring for air-leak and
intra-pleural pressure

- Patient data always visible in numeric and graphic mode

- Mini USB-key port for data downloading


Drentech Palm Evo

Kit portable chest drainage vacuum unit w/display _____ code 10192
Packaging: 1 unit for box

Chest unit protection filter _____ code 10538
Packaging: 36 units for box

Drentech Simple3 Evo

One patient tube _____ code 10125
Single use, sterile - Packaging: 4 units for box

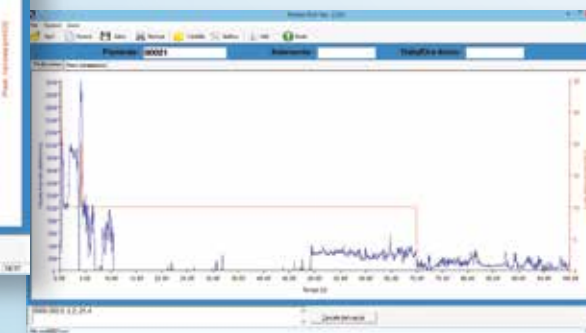

Improving drainage management provides the elements for a complete assessment of patient's status

- Total patient awareness: real time air leak data recording 99 hrs
- Real-time and graphic intrathoracic monitoring
- Immediately available in digital or graphic mode via touch screen mode
- Recording air leakage also in total gravity mode
- Equipped with video and audio allarms



Pictures kindly offered by: A.O.U. Città della Salute e della Scienza (TO)
Presidio Ospedaliero Molinette (Dr. P.L. Filosso)

DRENTECH™ mobile REDLINE

Mobile Redline VACUUM UNIT

TECHNICAL SPECIFICATIONS

- Negative pressure regulated through key-board from 0 up to -60 cmH₂O (0 to -6kPa)
- Current suction level always visible through display
- Autonomy of 36h

FEATURES and BENEFITS

- A complete freedom for the patient and an easy transport in hospital facilities
- Early patient mobilization and less work-load for the hospital staff
- Stable base to prevent accidental fall down
- High comfort for the patient (noiseless)

DISPOSABLE CHEST DRAIN

- Large capacity drainage unit
- Underwater seal valve separated from collection chamber
- All time clear visibility for air-leak directly on frontal
- Suitable for adult and paediatric use



DrenTech Mobile Redline

Kit portable chest drainage vacuum unit
Packaging: 1 unit for box

code 10188

Chest unit protection filter

Packaging: 36 units for box

code 10539



REDAX - Via Galileo Galilei, 18
46025 Poggio Rusco (MN) ITALY
Tel. +39 0386 830582 Fax +39 0386 51898
E-mail: info@redax.it - www.redax.it



REV.05-2019

DRENTECH™ REVOLUTION PROJECT

DRENTECH™ mobile REDLINE

DRENTECH™ PALMEVO

Independent

Mobile

Air-leak & pressure
monitoring

for Cardio-Thoracic
Surgery

